


**OLYMPUS**

Endo**Therapy**

Single Use Aspiration Needle

## **ViziShot 2**

Improved access and enhanced lymph node penetration  
for more efficient sampling in EBUS-TBNA

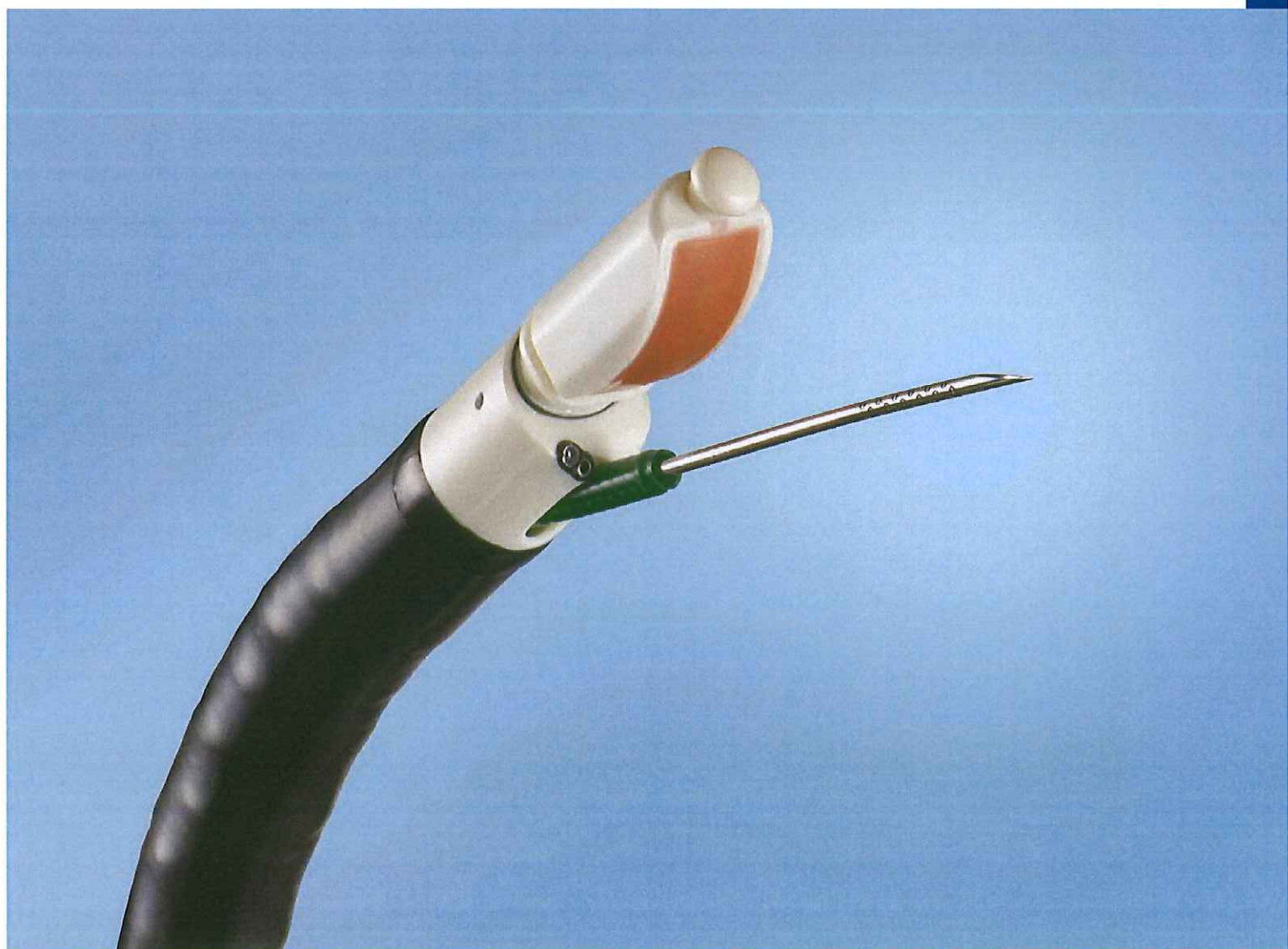


# SINGLE USE ASPIRATION NEEDLE

# ViziShot 2

## Features


### Improved flexibility

Advanced needle design allows improved flexibility, helping to support access to targets in bronchi which are usually difficult to reach with conventional needles.

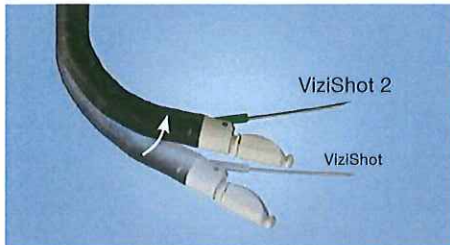


### Safety mechanism

ViziShot 2 incorporates a double-safety mechanism, the needle adjuster lever lock and the needle retract lock.

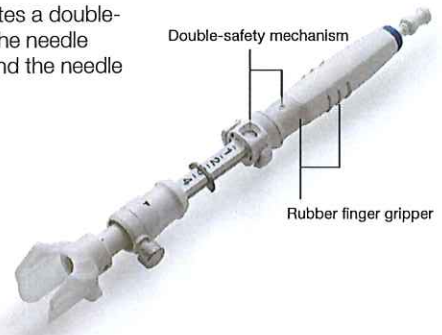


### Enhanced puncture capability

The sharper needle tip design reduces the puncture force required to take a sample.

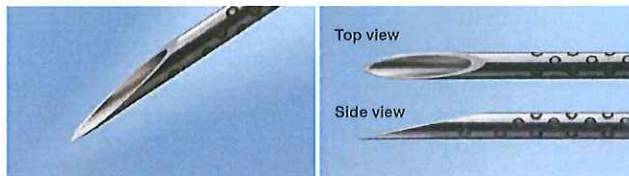

### Ergonomic design

New, easy-to-handle design with rubber finger gripper supports needle control and improves stability.

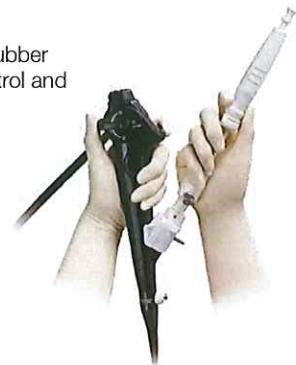

### Increased sheath visibility

Colored sheath tip increases visibility during endoscopic observation.




### Efficient echogenicity

Dimpled needle surface allows excellent ultrasound visibility.

## Specifications

Model Name	NA-U401SX-4021	NA-U401SX-4022
Maximum Insertion Portion Diameter	1.9 mm	
Working Length	700 mm	
Needle Gauge	21G	22G
Needle Length	Normal	20 mm
	Maximum	40 mm
Compatible OLYMPUS Ultrasonic Endoscope	Model and working length	BF (600mm)
	Channel inner diameter	2.0 mm/2.2 mm
	Ultrasonic scanning direction	Parallel to insertion direction

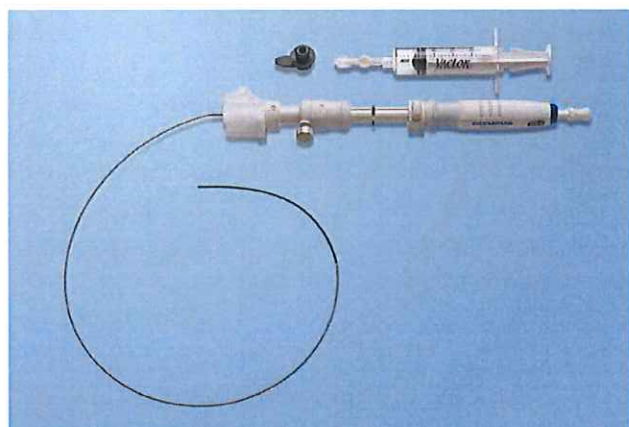

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

**OLYMPUS**

OLYMPUS MEDICAL SYSTEMS CORP.  
Shinjuku Monolith, 2-3-1 Nishi-Shinjuku, Shinjuku-ku, Tokyo 163-0914, Japan

For a complete listing of sales and distribution locations visit:  
[www.olympus.com](http://www.olympus.com)

Printed in Japan N8600387-042019