

Dilation Made Fast & Easy

- Rapid inflation, insertion and deflation times promote overall procedural speed and efficiency
- The low-profile balloon catheter is designed for quick and easy insertion and retrieval through the scope channel

Indications

- EZDilate Fixed Wire balloon is designed to dilate strictures in the Esophagus
- EZDilate Wire Guided Balloon is designed to dilate strictures in the Esophagus, Pylorus, Colon and Ampulla (Sphincter of Oddi with or without prior sphincterotomy, adults only)
- Pre-loaded stainless steel .035" guide wire supports dilation in tortuous anatomy
- The EZDilate Balloon Dilator is compatible with any 60 cc inflation device

ORDERING INFORMATION

EZDILATE FIXED WIRE BALLOON DILATORS, SINGLE USE – ESOPHAGEAL (1/BOX)

Part Number	Inflated Outer Ø (mm)	Inflated Outer Ø (Fr)	Inflation Pressure (atm)	Min. working channel Ø (mm)	Working length (mm)	Balloon length (mm)
BD-400P-0880	6.0 - 7.0 - 8.0	18-21-24	2.0-5.5-9.5	2.8	1900	90
BD-400P-1080	8.5 - 9.5 - 10.5	25.5-28.5-31.5	2.0-5.5-9.5	2.8	1900	90
BD-400P-1380	11.0 - 12.0 - 13.0	33-36-39	2.0-3.5-5.5	2.8	1900	80
BD-400P-1580	13.5 - 14.5 - 15.5	40.5-43.5-46.5	2.5-4.0-6.0	2.8	1900	80
BD-400P-1880	16.0 - 17.0 - 18.0	48-51-54	2.5-4.0-5.5	2.8	1900	80
BD-400P-2080	18.0 - 19.0 - 20.0	54-57-60	2.5-3.5-5.0	2.8	1900	80

EZDILATE WIRE GUIDED BALLOON DILATORS, SINGLE USE – ESOPHAGEAL, PYLORIC, COLONIC, BILIARY (1/BOX)


Part Number	Inflated Outer Ø (mm)	Inflated Outer Ø (Fr)	Inflation Pressure (atm)	Min. working channel Ø (mm)	Working length (mm)	Balloon length (mm)	Guidewire Ø (Preloaded & Compatible)
BD-410X-0855	6.0 - 7.0 - 8.0	18-21-24	2.0-5.5-9.5	2.8 mm	2400	65	.035"
BD-410X-1055	8.5 - 9.5 - 10.5	25.5-28.5-31.5	2.0-5.5-9.5	2.8 mm	2400	65	.035"
BD-410X-1355	11.0 - 12.0 - 13.0	33-36-39	2.0-3.5-5.5	2.8 mm	2400	55	.035"
BD-410X-1555	13.5 - 14.5 - 15.5	40.5-43.5-46.5	2.5-4.0-6.0	2.8 mm	2400	55	.035"
BD-410X-1855	16.0 - 17.0 - 18.0	48-51-54	2.5-4.0-5.5	2.8 mm	2400	55	.035"
BD-410X-2055*	18.0 - 19.0 - 20.0	54-57-60	2.5-3.5-5.0	2.8 mm	2400	55	.035"

* indicated for Esophageal, Colonic and Biliary only

INFLATION DEVICE FOR ENDOSCOPIC BALLOON DILATION, SINGLE USE (1/BOX)**

Part Number	Syringe Volume (cc)	Max Pressure (atm)
MAJ-1740	60	15

** Optional Accessory – EZDilate Fixed Wire and Wire Guided Balloon dilators are compatible with any 60 cc inflation device.

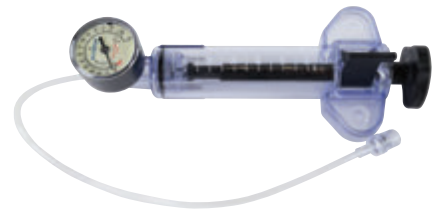

© 2015 Olympus America Inc. All rights reserved. Printed in USA. ™ & ® Trademark or Registered Trademark of Olympus or its affiliated entities in the U.S. and/or other countries of the world. Subject to change without notice. OAIURO0815BROxxxx

OLYMPUS AMERICA INC.
3500 Corporate Parkway, PO Box 610, Center Valley, PA 18034

For more information or order inquiries,
please call 800-848-9024 or visit us at
www.medical.olympusamerica.com

OLYMPUS®

Your Vision, Our Future

Endo**Therapy**

EZDilate

MULTI-STAGE ENDOSCOPIC BALLOON DILATOR

Controlled Accurate Dilation: Reliability in Exact Sizing

Precise Placement: Confidence in Positioning with Optimal Clarity

Versatility in Design & Performance: Expand your procedural options



CONTROLLED RADIAL DILATION WITH GREATER PRECISION

EZDILATE THREE STAGE BALLOON DILATORS OFFER GRADUAL AND ACCURATE DILATION

Indicated for strictures throughout the GI tract—including dilation of the Sphincter of Oddi (Endoscopic Papillary Balloon Dilatation/EPBD).

EZDilate facilitates precise placement in the stricture due to its unique, patent pending center marker as well as excellent visualization through the balloon during dilation for improved monitoring of the stricture.


FIXED WIRE: indicated for dilation in the Esophagus



WIRE GUIDED: indicated for dilation in the Esophagus, Pylorus, Colon & Papilla



Offered in both FIXED WIRE and WIRE GUIDED designs, with sizes ranging from 6-7-8 mm to 18-19-20 mm

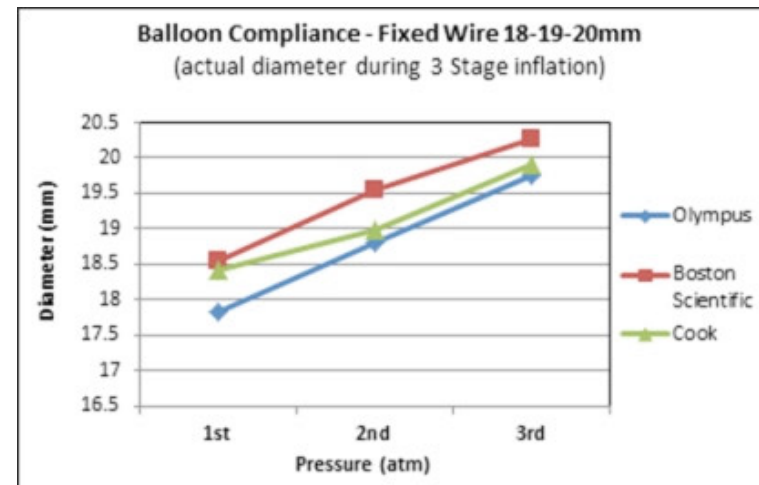

Controlled Accurate Dilation


- Engineered to provide superior diameter accuracy and precision, without the need to over inflate to reach the target specification
- Robust nylon material disperses radial pressure along the full length of the balloon ensuring repeatable, reliable balloon dilation with minimal waisting
- Semi compliant balloon designed to maintain diameter and shape in all stricture types

Versatility in Design and Performance

- Portfolio range to suit your procedure needs throughout GI Tract—including Biliary application
- Unique and accurate sizing provides the opportunity to inflate at half millimeter intervals
- Balloons Compatible with the 60 cc inflator of choice—providing flexibility in clinician preference and overall ease of use

Precise Placement: Confidence in Positioning with Optimal Clarity

- Low Profile catheter design and specialized balloon coating promotes smooth and efficient insertion with minimal resistance
- Patent Pending Center Visual Marker facilitates placement of balloon, positioning accuracy & assessment of stricture during dilation
- Radiopaque markers provide enhanced visualization under fluoroscopy
- Ultra clear nylon material enables enhanced visualization of the stricture through the balloon during dilation

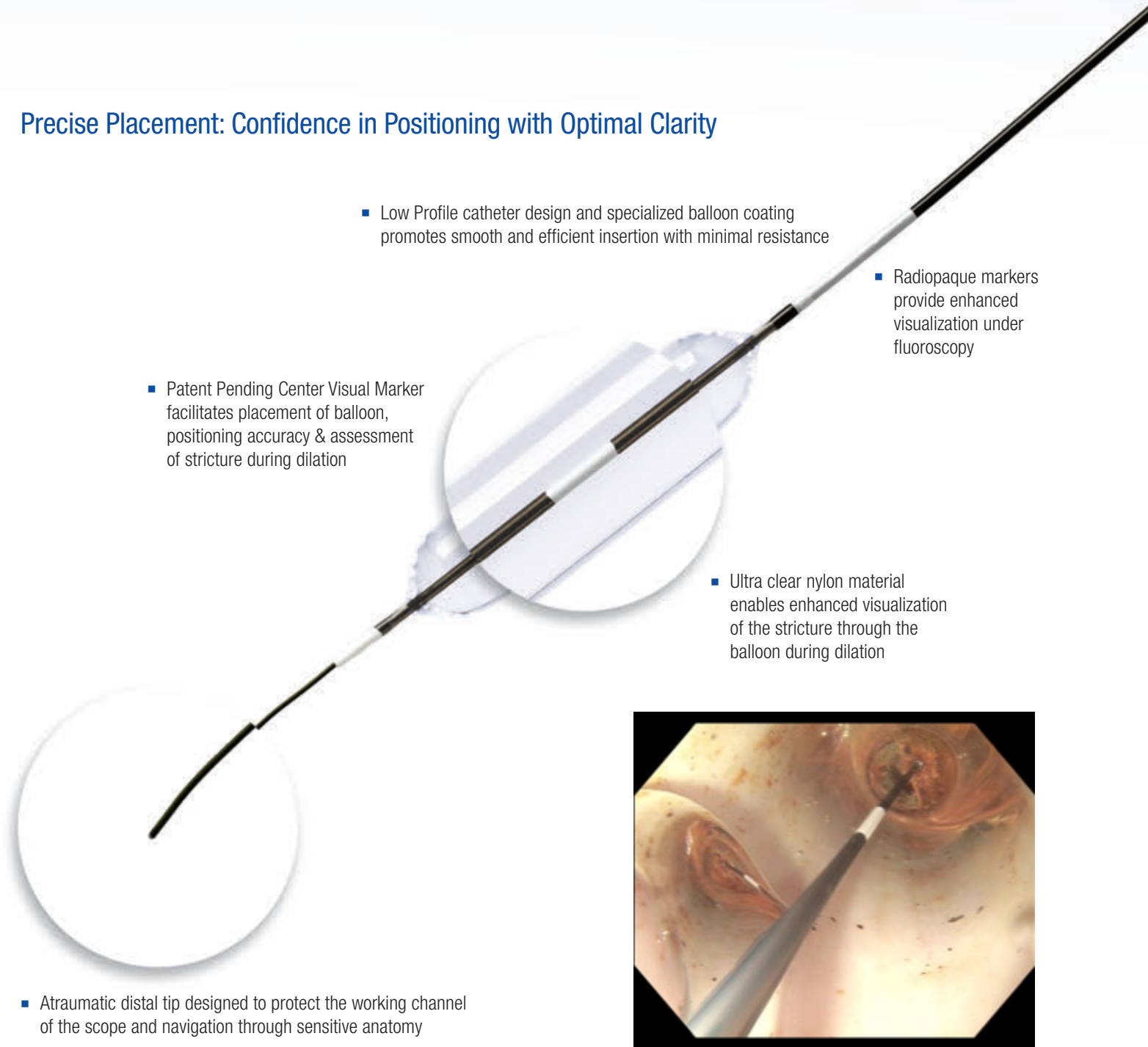

- Atraumatic distal tip designed to protect the working channel of the scope and navigation through sensitive anatomy